

# Accessories

for Merivaara Operating Tables



# Accessories

## Comprehensive range of operating table accessories

This comprehensive selection of accessories is designed to complement Merivaara operating tables in a variety of surgical procedures. Operating table accessories represented in this brochure are manufactured or distributed by Merivaara. When equipped with these high-quality accessories the tables are also suitable for many specialized procedures. The safe and easy fixing and ergonomic design ensure a smooth surgical process with less time demands. All materials used have been carefully selected. The soft, body conforming surfaces provide comfort and safety for the patient, as well as meeting surgical requirements.



## Contents

Arm rests .....	4
Arm and hand tables .....	4
Supports and straps.....	5
Head rests.....	6
Leg rests .....	7
Leg sections.....	7
Leg holders and supports.....	8
Rail extensions .....	9
Anesthesia screens and infusion rod.....	9
Clamps .....	10
Bowls.....	10
Control devices .....	11
Table top sections.....	11
Carbon fibre sections.....	11
Table top sections for Promerix.....	12
Table top sections for Practico.....	12
Pads.....	14
Other accessories .....	15
Special devices.....	15
Orthopedic extension device .....	16
Recommended accessories for the most common fields of surgery .....	18

# Accessories

## Arm rests



### 1. Arm rest 19180

Padded arm rest, incl. clamp and strap. Adjustable height and angle (ball joint), easy to operate and adjust by one lever.

Order no. 1000 19180  
1000 19180U for US-rail

Suitable rails: EUR, US



### 2. Arm rest with radial clamp and ball joint 19185

Padded arm rest, incl. clamp and strap. Adjustable height and angle (ball joint), easy to operate and adjust by one lever. Fixing to rail with radial clamp to give wider adjustment angles. Incl. radial clamp and strap 110.

Order no. 1000 19185  
1000 19185U for US-rail

Suitable rails: EUR, US



### 3. Arm rest 125

Padded, includes a fixing clamp and straps 110. Adjustable height. Gas spring assisted angle adjustment.

Order no. 1000 00125

Suitable rails: EUR, US



### 4. Easy armboard 11801

Easy to install to the siderail. Spring loaded design simply snaps onto the rail and is released by a convenient lever accessible from both sides. Light weight armboard rotates 180°. Fully radiolucent with comfortable pad and velcro strap.

Order no. 1000 11801(EUR)  
1000 11801B (UK)

Suitable rails: EUR, UK



### 5. Lateral armboard 11802

Provides arm support for prone, neuro, lateral and chair positions. A lockable ball joint combined with 150 mm of horizontal travel permits a wide range of position flexibility. The premium durasoft pad ensures patient comfort. Universal clamp 10300 is needed for fixing.

Order no. 1000 11802 Lateral arm board  
1000 10300 Universal clamp

Suitable rails: EUR, US, UK



### 6. Arm and hand table 11911

Radiolucent top (350 x 800 mm), light to handle. Swivel-ball leg slides, angles and folds provide unimpeded access. Adjust vertically for tables with X-ray tops. Endrail provided for traction attachment. Mounts to siderail with two included clamps 11043 (EUR/UK). For US-rail 2 pcs clamps needs to be ordered separately. Soft, premium pressure management pad improves patient comfort and stability.

Order no. 1000 11911 Arm table  
1000 11043 U Clamp for US-rail

Suitable rails: EUR, US, UK

## Arm and hand tables



### 7. Carpa table 11912

A reduced-size (350 x 650 mm) arm and hand table for carpal tunnel and other minor procedures. Light weight, 2.7 kg. Unique cantilevered side frames create an extremely rigid surface that is completely radiolucent. Vertically adjustable 50 mm. Soft pressure management pad for optimum patient comfort and stability. Mounts to siderail with two included clamps, 11043 (EUR/UK). For US rail 2 pcs clamps needs to be ordered separately.

Order no. 1000 11912 Carpa table  
1000 11043U Clamp for US-rail

Suitable rails: EUR, US, UK



### 8. Arm operating table 10390

Special shaped arm operating table with one support leg. Option to use arm traction device 10392 and hand fixation device 10393 ("Weinberger - type"). Includes mattress and clamps.

Order no. 1000 10390 Hand operating table 10390  
1000 10392 Traction device 10392 for 10390 only  
1000 10393 Handfixation device 10393 for 10390 only  
1000 10390U Hand operating table 10390U for US-rail

Suitable rails: EUR, US



### 9. The Weightless Shoulder Suspension System 64436

Includes one sterile Arm Trap Disposable Sterile Arm Trap 1000 21200 for shoulder suspension system, arm trap is affixed to patient with CoFlex® NL Self-Adherent Wrap (roll included). Incl. 6 pcs. Clamp 10300 is required for fixing.

Order no. 1000 64436 Weightless shoulder suspension system  
1000 10300 Universal clamp  
1000 21200 Disposable Sterile Arm Trap

Suitable rails: EUR, US, UK



### 10. Humerus traction device

Traction device 11443 for intramedullary nailing of the humerus, (requires clamp 10300). Padded counteraction post for humerus, (requires clamp 10300).

Weinberger hand traction device, for hand and finger traction.

Order no. 1000 11443 Traction device 11443  
1000 11438 Counteraction post  
1000 10300 Clamp 10300  
1000 10393 Weinberger hand traction device

Suitable rails: EUR, US, UK





# Supports and straps



**11. Wrist straps (pair) 107**  
 Fixed to the accessory rail with clamps. Safety belt type.  
 Order no. 1000 00107  
 Suitable rails: EUR



**12. Leg / Body strap 10886, wide**  
 Fixing with velcro-tape.  
 Order no. 1000 10886  
 Suitable rails: EUR, US, UK



**13. Patient restraint strap 11081**  
 Attaches to the siderail quickly with stainless steel hooks. Adjust easily with two airplane-style buckles. The strap is made of strong, conductive, reinforced rubber for durability and ease of cleaning.  
 Order no. 1000 11081  
 Suitable rails: EUR, US, UK



**14. Side support 120**  
 Includes a clamp and PU-pad. Adjustable height, angle and distance to patient.  
 Order no. 1000 00120  
 Suitable rails: EUR



**15. Lateral braces 11651 and 11652**  
 11651 with round pad, 11652 with rectangular pad. Provide secure and comfortable positioning for hip and other procedures requiring optimum lateral stability. The durossoft pads provide soft cushioning of bony prominences. The swivelling heads conform to anatomical contours, can be securely locked for precise positioning. Includes mounting clamps 11043. For US-rail 2 pcs clamps needs to be ordered separately.  
 Order no. 1000 11651 Lateral brace, round pad  
 1000 11652 Lateral brace, rectangular pad  
 1000 11043 U Clamp for US-rail  
 Suitable rails: EUR, US, UK



**16. Supporting Roll, vertical**  
 Order no. 1000 10368 (requires clamp 10301)  
 Suitable rails: EUR



**17. Supporting Roll, horizontal**  
 Order no. 1000 10369 (requires clamp 10301)  
 Suitable rails: EUR



**18. Shoulder supports 19129 (pair)**  
 Includes clamps. Padded. Adjustable angle and distance to shoulders.  
 Order no. 1000 19129  
 Suitable rails: EUR



**19. Shoulder braces 11295 (pair)**  
 Provides comfortable support for Trendelenburg positioning. The durossoft pads on broad support surfaces tilt as needed to conform to bony shoulder prominences and minimize pressure. For fixing Rail clamp 11043 is needed, 2 pcs. (not included).  
 Order no. 1000 11295  
 Suitable rails: EUR (11043), US (11043U), UK (11043)



**20. Straps for armrest 110 (pair)**  
 Straps for 110 for arm rests 125, 19180 and 19185.  
 Order no. 1000 00110




**21. Support belt**  
 Support belt for Beach Chair.  
 Order no. 1000 20127  
 Incl. a clamp, 1 pcs  
 Suitable rails: EUR, US, UK



**22. Leg- / arm protector 60485**  
 Padded leg- / arm shield for safe patient positioning.  
 Order no. 1000 60485


# Accessories

## Head rests

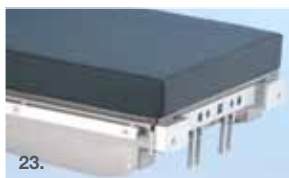

**23. Fixing adapter 60146**  
For Doro/Mayfield® Cranial Positioning System and other special head rests. For Promerix, Practico.  
Order no. 1000 60146

Suitable with: [Promerix](#), [Practico](#)

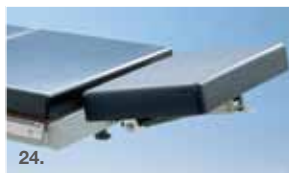

**24. Wide headrest 141**  
Size 270 × 400 mm. Fits all Merivaara operating tables. For Promerix and Practico the fixing adapter is needed.  
Order no. 1000 00141

Suitable with: [Promerix](#), [Practico](#), [Rapido](#), [OP1650](#)



**25. Narrow headrest 142**  
Size 200 × 270 mm. Fits all Merivaara operating tables. For Promerix and Practico the fixing adapter is needed.  
Order no. 1000 00142

Suitable with: [Promerix](#), [Practico](#), [Rapido](#), [OP1650](#)



**26. Head rest 18101 (wide)**  
With Rapido for general surgery, all Merivaara operating tables with "Beach Chair"- back rest. For Promerix and Practico the fixing adapter is needed.  
Order no. 1000 18101.

Suitable with: [Promerix](#), [Practico](#), [Rapido](#), [OP1650](#)



**27. Special head rest 18143**  
For ENT-surgery. Fits all Merivaara operating tables. For Promerix and Practico the fixing adapter is needed.  
Order no. 1000 18143

Suitable with: [Promerix](#), [Practico](#), [Rapido](#), [OP1650/1700](#)



**28. Special head rest 18150 for ophthalmology & ENT**  
For ophthalmology and ENT-surgery. Fits to all Merivaara operating tables. For Promerix and Practico the fixing adapter is needed.  
Order no. 1000 18150

Suitable with: [Promerix](#), [Practico](#), [Rapido](#), [OP1650/1700](#)



**29. Special head rest 18151 for ophthalmology**  
For Ophthalmology, double joint mechanism, wider range of adjustments. For Promerix and Practico the fixing adapter is needed.  
Order no. 1000 18151

Suitable with: [Promerix](#), [Practico](#), [Rapido](#), [OP1650/1700](#)



**30. Special head rest 30126**

For "Beach Chair" back rest for shoulder and upper arm procedures. Unique sliding function with double joint mechanism, very wide adjustments in all directions.

Fits to Promerix 60525 and Practico 30125 back rests.

Order no. 1000 30126.

Suitable with: [Promerix](#), [Practico](#), with [Beach chair back rests](#)



**31. Helmet type head rest 30128**

Helmet type head rest 30128 for shoulder and upper arm procedures with "Beach Chair" back rest.

Unique sliding function with double joint mechanism, very wide adjustments in all directions. Fits Promerix 60525 and Practico 30125 back rests.

Order no. 1000 30128

Suitable with: [Promerix](#), [Practico](#), with [Beach chair back rests](#)



**32. Neurosurgical head rest 144**

Padded. Fits all Merivaara operating tables. For Promerix and Practico the fixing adapter is needed.

Order no. 1000 00144

Suitable with: [Promerix](#), [Practico](#), [Rapido](#), [OP1650](#)



**33. Arm supports 151 (pair)**

Fixing to the accessory rails. Includes clamps.

Order no. 1000 00151

Suitable rails: [EUR](#), [UK](#), [US](#)



**34. Push bar to the head end of Promerix and Practico 20593**

Suitable for 60140 and 60142 Head rests.

Order no. 1000 20593

Suitable with: [Promerix](#), [Practico](#)



**35. Gas spring assisted head rest 60140**

Optional mattresses moulded IS or visco elastic VEF. Accessory rail option EUR or US.

Order no. 1000 60140

Suitable with: [Promerix](#), [Practico](#)



**36. Double articulation head rest 60142**

Gas spring assisted. Optional mattresses moulded IS or visco elastic VEF. Including accessory rails (EUR, US).

Order no. 1000 60142

Suitable with: [Promerix](#), [Practico](#)

## Leg rests



**37. Leg rests 20118 (pair)**  
Padded and with velcro straps and a ball joint. Stainless steel. Universal clamp 10300 needed for fixing.

Order no. 1000 20118

10300 for all rails

37.



**38. Powered Leg rests 60482 (pair) for Promerix table's powered leg joint**

Provides excellent patient comfort and precise surgical positioning. The VEF-pad reduces popliteal pressure. Velcro straps for security. Integrated clamps included.

Order no. 1000 60482

Suitable with: Promerix

38.



**39. Knee crutches 11181 (pair)**

Provides excellent patient comfort and precise surgical positioning. The "Doughy soft" pad reduces popliteal pressure. Velcro straps for security. Universal clamp 10300 is needed for fixing.

Order no. 1000 11181

10300 for all rails

39.




**40. Stirrup system 12581 (pair)**

Eliminates pressure of the popliteal fossa and perineal nerve. Eliminates uncomfortable rotation, stretching and compression of the knee and hip joints. Universal clamp 10300 is needed for fixing.

Order no. 1000 12581

10300 for all rails

40.

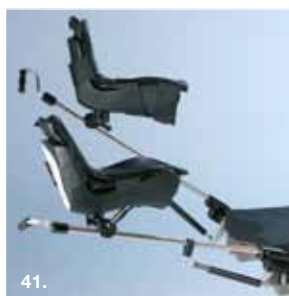

**41. Gas spring assisted, stirrups 12583 (pair)**


Allows direct placement with siderail clamps and light weight molded boots which are lined with a durable, soft pad for optimum patient comfort and protection. Each stirrup joint has lithotomy degree indicators. Clamps 11043 are needed for fixing.

Order no. 1000 12583 Stirrups, pair  
1000 11043 Clamps, 2 pcs.  
1000 11043 U Clamps for US-rail, 2 pcs.

Suitable rails: EUR (11043), US (11043U), UK (11043)

41.

## Leg sections

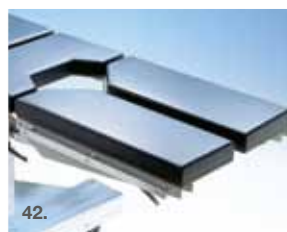

**42. Divided leg section**

Section replaces one-piece foot section. Gas spring assisted.


Order no. 1000 00148 OP1650

1000 18261 Rapido

1000 30248 Practico, Promerix

Suitable with: Rapido, Practico, Promerix, OP1650

42.

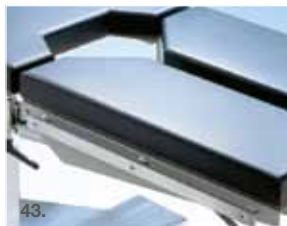

**43. Accessory rails for divided leg sections**

Optional rail for Rapido/ 1650 divided leg sections.

Order no. 1000 01481 1650/Rapido, for EUR-rail

1000 01482 Rapido, UK-rail

1000 01483 Rapido, US-rail

Suitable with: Rapido, OP1650

43.



**44. Foot rest/table extension 20114**

Suitable for all Merivaara tables. Available both IS and VEF mattress, or without the mattress. Universal clamps (2 pcs) included.

Order no. 1000 20114

Suitable with: all Merivaara tables

44.

# Accessories

## Leg holders and supports



### 45. Leg holder 12551 (pair)

Stainless steel construction for smoother operation and improved reliability. Strong rotation lock prevents slippage and permits vertical adjustment from 700–1030 mm. Universal clamp 10300 is needed for fixing.

Order no. 1000 12551

Suitable rails: EUR, UK, US



### 50. Arthroscopic stress post 12530

A levering fulcrum for improved access to the medical compartment of the knee. Mounts to accessory rail. Universal clamp 10300 is needed for fixing.

Order no. 1000 12530

10300 for all rails



### 46. Thigh holder 12012

Rotating, complete with antistatic pad. Universal clamp 10200 is needed for fixing.

Order no. 1000 12012 (incl. 1 pc)

10200 for all rails



### 51. Total knee stabilizer 12585

Convenient, stable device to support the flexed leg during total knee procedures. Includes pad and clamp 11043 (EUR/UK) For US-rail clamp needs to be ordered separately. Extra pads 12586 are available.

Order no. 1000 12585 Total knee stabilizer

1000 12586 Pads, extra pieces

1000 11043U Clamp for US-rail

Suitable rails: EUR, UK, US



### 47. The prepper 12552

The padded cradle comfortably supports patients leg during prep procedure. Shaft can be rotated to move offset in various positions for convenient access to the leg. Generous 750 mm long, 16 mm post. Universal clamp 10300 is needed for fixing.

Order no. 1000 12552

10300 for all rails



### 52. Arthroscopic Leg holder 10353

Needs clamp 1000 10308 for fixing.

Order no. 1000 10353

1000 10308 Clamp

Suitable rails: EUR, US, UK



### 48. Light foot rest

Provides a light-weight, easy-to-use foot rest, that supports patient's legs prior to gynaecological or urological procedures (transport, anaesthesia).

Order no. 1000 18001 for Rapido,

1000 20201 for Practico and  
Promerix

1000 12001 for OP1650

Suitable with: Promerix, Practico, Rapido,  
OP1650



### 49. Light foot rest

Provides a light-weight, easy-to-use foot rest with latch locking mechanism. Supports patient's legs prior to gynaecological or urological procedures (transport, anaesthesia).

Order no. 1000 60201

Suitable with: Promerix, Practico



## Rail extensions




### 53. Accessory rail extension 19128

Stainless steel, length 550 mm. Fits to both sides of the table. Suitable for all Merivaara operating tables with European standard siderails.

Order no. 1000 19128  
1000 19128U

Suitable rails: [EUR \(19128\)](#), [US \(19129U\)](#)

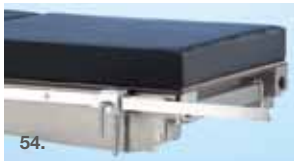


### 54. Accessory rail extension 10355

Length 400 mm

Order no. 1000 10355

Suitable rails: [EUR \(19128\)](#), [US \(19129U\)](#)

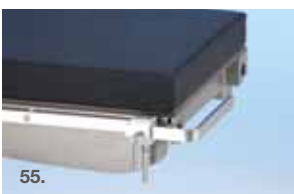

### 55. Accessory extension rail, angled 10356

Angeled, Length 270/150 mm

Order no. 1000 10356

Suitable rails: [EUR \(19128\)](#), [US \(19129U\)](#)



### 56. Angular accessory rail 10359

Length 250 mm

Order no. 1000 10359

Suitable rails: [EUR \(19128\)](#), [US \(19129U\)](#)

## Anesthesia screens and infusion rod



### 57. Anaesthesia screen 111

Adjustable height and angle. Made of stainless steel. Universal clamp 10300 is needed for fixing.

Order no. 1000 00111

### Extension arm 116 (single)

An extra wing is attached to the anaesthesia screen 111, one or two pcs.

Order no. 1000 00116

[10300 for all rails or 10308](#)



### 58. Flexible anaesthesia screen 11001

Malleable screen conforms to any shape desired. Universal clamps 10200 (2 pcs) are needed for fixing.

Order no. 1000 11001

[10200 for all rails](#)



### 59. Infusion rod 749

Chrome plated, fixed to the accessory rail. Includes fixing clamp. Adjustable height, two or four hooks.

Order no. 1000 00749

Suitable rails: [EUR](#), [US](#), [UK](#)



### 60. Infusion rod holder 18578

Fixed to the base casing, storage holder for Rapido, infusion rod with two hooks.

Order no. 1000 18578

Suitable rails: [Rapido](#)

# Accessories

## Clamps



### 61. Accessory Clamp 10200 Universal Ø 16

Accessory Clamp, Universal, ø 16  
Universal model, fits all railtypes.  
Aluminium.

Catch the ø 16 round pole in vertical or in horizontal direction and Eur-rail size flat bar in horizontal direction.

Order no. 1000 10200

Suitable rails: EUR, US, UK



### 62. Accessory Clamp 10301 Ø 20 mm

Accessory clamp 10301, Ø 20 mm, for EUR-rails.

Order no. 1000 10301

Suitable rails: EUR



### 63. Radial Accessory Clamp 10300 Universal Ø 16-20

Radial Accessory Clamp ø16-20, Universal model, fits all railtypes. Stainless steel. Easy to pre-adjust for different diameters ø16-20. Pat. Pend. Extremely good grip to round pole, particularly handy with the leg rests and stirrups with round pole.

Order no. 1000 10300

Suitable rails: EUR, US, UK



### 64. Rail clamp 11043

Secure locks all accessories with flat mounting posts anywhere along the table siderail. Perfect for table top extensions and knee board applications. Stronger, lighter weight replacements for original rail locks. Reduces the free play associated with standard clamps.

Order no. 1000 11043

1000 11043U

Suitable rails: EUR, UK (11043), US (11043U)



### 65. Radial Clamp 10308

For round pole, up to Ø 20 mm. Fits to all rails EUR, US, UK (not for leg rests and stirrup systems).

Order no. 1000 10308

Suitable rails: EUR, US, UK

## Bowls



### 66. Bowl 6 I

201571 without drainage, 201572 with drainage. Stainless steel (6 l). Fixed to siderail with a clamp. Fits all Merivaara operating tables.

Order no. 100 201571 or 100 201572

### Bowl 6 I, for Promerix, Practico and Rapido without attachment

Bowl without drainage 20501, bowl with drainage 20502

Order no. 1002 01501

1002 01502

Suitable rails: EUR, US



### 67. Bowl 6 I, for Promerix

61581 without drainage, 61582 with drainage. Stainless steel. Fixed under the table top. Holder for 6l bowl, fixed under the base, Promerix 61580.


Attached to gyn-section 60440 or 60425.

Order no. 1000 61580 Holder

1000 61581 Bowl without drainage + holder

1000 61582 Bowl with drainage + holder

Suitable with: Promerix

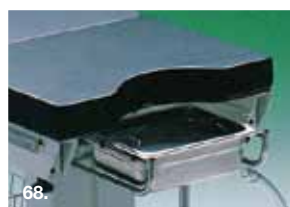

### 68. Bowl 6 I, for Rapido

181581 without drainage, 181582 with drainage. Stainless steel. Fixed under the table top.

Order no. 1001 81581 bowl without drainage

1001 81582 bowl with drainage

Suitable with: Rapido



### 69. Bowl 4 I, for 1650

158 without drainage. Stainless steel. Fixed under the table top.

Order no. 1000 00158

Suitable with: OP1650



### 70. Screen 18159

Separately for bowls, 6 l.

Order no. 1000 18159

Suitable with: Bowls 6 I

## Control devices



71.

### 71. Control devices for Promerix

Hand control for Promerix

Order no. 1000 60800 Hand control, wired  
1000 60820 IR-hand control

Suitable with: Promerix




72.

### 72. Merimote

Merimote™ is a multifunctional remote control for Merivaara's Promerix operating table. With intuitive touch screen Merimote™ guides to get full advantage of the working environment and gives access to more advanced product support.

Order no. A42682400

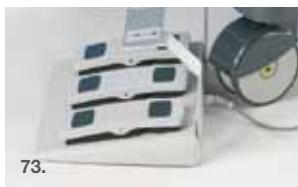

73.

### 73. Foot control for Promerix & Practico

Foot control device Height, Trend., Tilt (optional Back on sliding).

Order no. 1000 60850

Suitable with: Promerix, Practico

## Table top sections



74.

### 74. Gynaecologic/urologic extension 18153 for Rapido

Gives more space (+150 mm) for the surgeon and makes C-arm use easier. Includes mattress.

Order no. 1000 18153

Suitable with Rapido



75.


### 75. Table width extender 10659, 10659U (single)

Increases the support surface for large patients. Can be used individually or in multiples to create a custom fit. Size 100 x 500 mm. Fixed to the accessory rail. Includes pad and accessory rail.

Order no. 1000 10659 (EUR)  
1000 10659U (US)

Suitable with: Promerix, Practico

## Carbon fibre sections

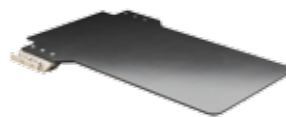

76.

### 76. Carbon Fibre Back-/Leg Plate 60500


Promerix Carbon Fibre section can be used as a back rest or as a leg plate. Length 950 mm (= like 40-section +25-section + head plate), replaces normal Promerix back and/or leg sections. Width 540 / 450 mm. Max. patient weight up to 275 kg. Length of 360-deg. imaging area 650 mm.

Mattress for the Carbon Fibre back-/leg plate 60501.80 mm VEF (Visco elastic foam), antistatic, welded seams.

Order no. 1000 60500 Carbon fibre back-/leg plate

1000 60501 mattress for carbon fibre back plate

Suitable with: Promerix

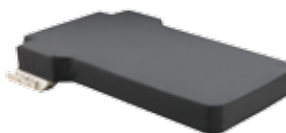

77.

### 77. Carbon Fibre Arm Rest 60510

Specially made to used with the Carbon Fibre back rest. The patient's counterweight is used to hold this arm rest on its palace. No metal parts, 100% Carbon Fibre – construction.

Total length 900 mm, arm rest size 160 x 680 mm.

Order no. 1000 60510

Suitable with: all tables



78.


### 78. Mattress for the Carbon Fibre Arm Rest 60511

Pad for the Carbon Fibre Arm Rest 80 mm VEF (Visco elastic foam), antistatic, welded seams. W = 160mm, L = 600 mm.

Order no. 1000 60511

# Accessories

## 79. a Table top sections for Promerix

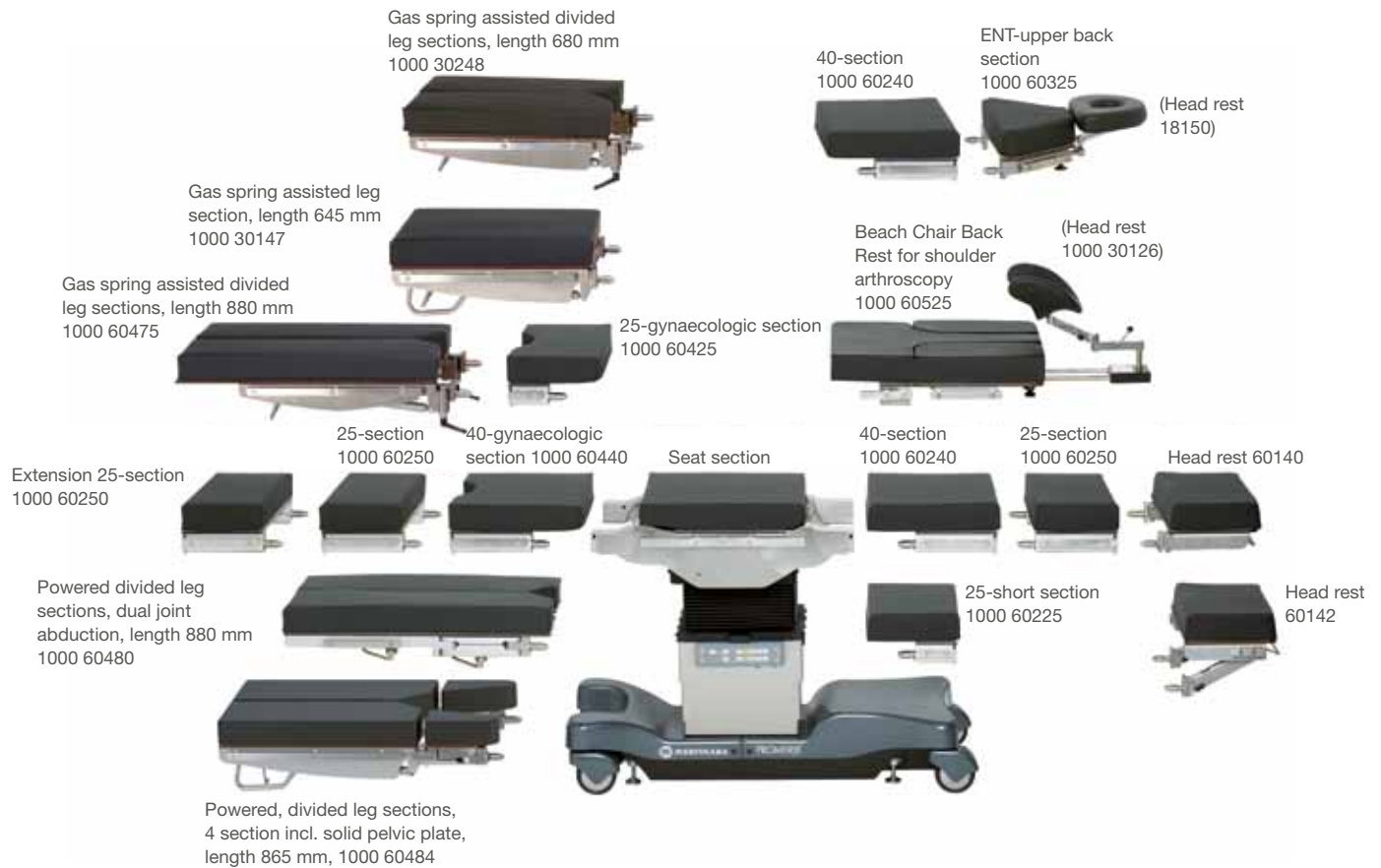

## 79. b Table top sections for Practico

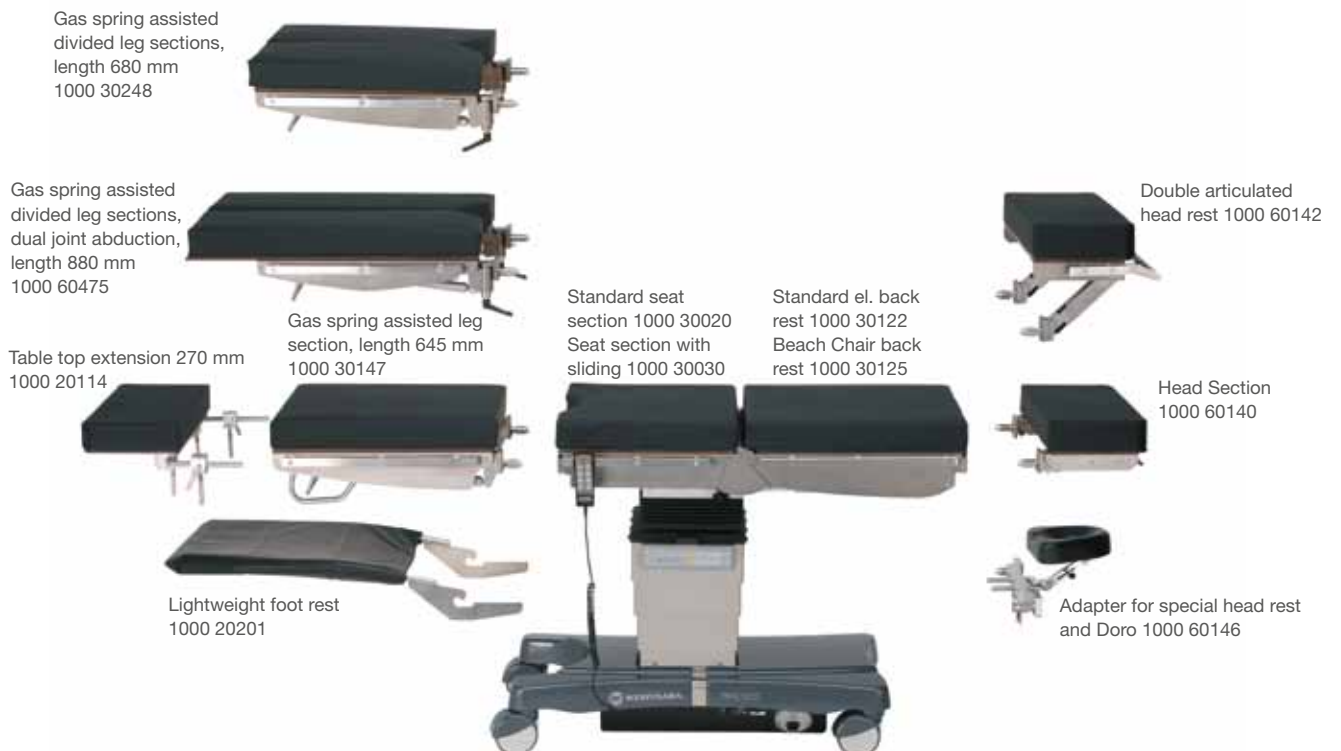




Table top section		Promerix: Section code	Practico: Section code	Mattresses VEF 80 mm	Mattresses IS 65 mm
Head rests	Gas spring assisted headrest	100060140	100060140	100060951	100060901
	Double articulation headrest	100060142	100060142	100060951	100060901
	Head Rest for Kidney and Thorax position	100060147	100060147	100060951	100060901

## BACK SECTIONS

Lower back rest	25-short section	100060225		100060955	100060909
	40-section	100060240		100060956	100060908
Upper back rest	25-section	100060250		100060955	100060902
	ENT-upper back section	100060325		100060962	100060912
	Beach Chair backrest	100060525		upper: 1000 60963 set	upper: Not available
				lower: 1000 60955	lower back: 1000 60909
Back sections: Practico	Beach Chair backrest, Practico, - factory assembled		100030125	upper: 1000 60963 set	Not available
				lower: 1000 60955	Not available
	Back rest, electrical, Practico, - factory assembled		100030122	A4220 9700	A4228 4100

## LEG SECTIONS

Upper leg section / Pelvic support	40-section	100060240		100060956	100060908	
	25-gynaecologic short section	100060425		100060959	100060907	
	40-gynaecologic section	100060440		100060957	100060905	
	Divided leg section, dual joint abduction		100060480		right: 1000 60961	right: 1000 60911
					left: 1000 60960	left: 1000 60910
		Divided leg plates, 4-sectional, incl. solid pelvic plate	100060484		Pelvic A4232 7400	Not available
				Femur A4235 0900, pair	Not available	
				Legs A4233 0100, pair	Not available	
Lower leg section	Divided leg section, gas spring assisted, long	100060475	100060475	1000 60961 right	1000 60911 right	
				1000 60960 left	1000 60910 left	
	Divided leg section, gas spring assisted, short	100030248	100030248	Leg L&R A4233010, pair	Leg L: A4230 5900 Leg R: A4230 5800	
	One piece leg section, gas spring assisted	100030147	100030147	A4220 9700	A4228 4100	
	25-section	100060250		100060955	100060902	
	Lower leg section, incl. 2x25-section	100060465		1000 60955, 2 pcs	1000 60902, 2 pcs	
Seat sections	Seat section without gynaecological cutout			100060953	100060903	
	Seat section with gynaecological cutout			100060954	100060904	
Seat + Back section	Seat section with gynaecological cutout + back section, for kidney elevator back section			1000 60952		

Available with moulded IS-mattress or welded visco-elastic memory foam mattress. Rail options: EUR, UK or USA.

## 80. PROMERIX VEF-mattress options separate/ Long seat-back mattress



25-gyn. + seat + 40-section



Seat with gyn. + 40-section + 25-section



40-gyn. + seat + 25-short section



Long seat + back mattress

# Accessories

## Pads

### 81. Stratus gel pads and positioners



**Closed head ring**  
size 140x60x35mm code 1000 66002  
size 200x75x45mm code 1000 66003



**Open head ring**  
size 140x60x35mm code 1000 66012  
size 200x75x45mm code 1000 66013



**Head pad**  
size 270x230x70 mm code 1000 66110



**Contoured headrest**  
size 190x190x60 mm code 1000 66081



**Prone headrest**  
size 280x230x145 mm code 1000 66111



**Heel support pads, pair**  
size 200x110x70 mm code 1000 66031



**Patient positioner**  
380x125x63 mm code 1000 66096



**Chest roll with flat bottom surface**  
405x150x100 mm code 1000 66042



**Armboard pad**  
600x125x20 mm code 1000 66072



**Crutch stirrup pads, 2 pcs**  
size 315x220x10 mm code 1000 66801



**Table pads**  
size 500x500x10 mm code 1000 66021  
size 1170x500x10 mm code 1000 66022



### 82. D-pillow 19497, visco-elastic ("memory") foam

Includes fixing straps, 21 x 49 x 10 cm.  
Order no. 1000 19497



### 85. Transfer mattress 19400

Used for transferring a patient between different patient platforms. Easy to use, easy to wipe clean, waterproof material. Colour dark blue.

Size 180 x 55 cm, CE-marked

Order no. 1000 19400



60571

83./84.

### 83. Prone Positioning System 60750

For spine and back procedures in prone position. 4 pcs of pads, including formed head support.

Width 540 mm

Recommended option + 84. D-pillow 19497

Order no. 1000 60750

### 84. Head Pillow for Prone position separated 60751

Order no. 1000 60751

## Other accessories




### 86. Push bars to foot end 18168, pair, foldable

Fits to all tables with EUR-rail in the leg section.

Order no. 1000 18168

Suitable rails: EUR

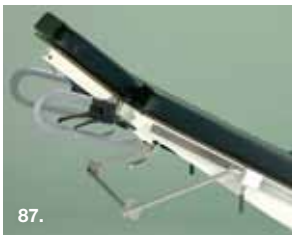

### 87. Paper roll holder for Rapido 18572

All rail types. Clamps included.

Order no. 1000 18572

Suitable rails: EUR, US, UK



### 88. Side rails 18173/18174

For all tables with EUR rail.

Order no. 1000 18173 chrome plated  
1000 18174 epoxy coated

Suitable rails: EUR



### 89. Wire Basket for utilities 18183

For Rapido.

Order no. 1000 18183

Suitable with: Rapido



### 90. Oxygen bottle holder 18189

Order no. 1000 18189 for 5l bottle

Suitable with: Rapido



### 91. Cassette tray

Suitable for different X-ray cassette sizes. Slides on the X-ray cassette rails under the table. Stainless steel.

Order no. 1000 18126 for Rapido  
1000 00126 for Promerix, Practico, OP1650

Suitable with: Promerix, Practico, Rapido, OP1650



### 92. Accessory stand 19220

For storing operating table accessories. Stainless steel. Fitted with accessory rails and steel baskets (4 pcs).

Order no. 1000 19220 or 1000 19220U (US rail)

Suitable rails: EUR, UK, US

## Special devices



### 93. Back and kidney elevator 19280

To be placed on the table top or on the mattress. The angle is adjustable with a crank mechanism. A separate mattress is included.

Order no. 1000 19280

Suitable with: Promerix, Practico, OP1650



### 94. Proctology accessory 154/30154

Allows the surgeon good access to the patient. SWL 180 kg

For Promerix 60425 (or 60440) needed.

Order no. 1000 30154 Practico & Promerix  
1000 00154 OP1650

Suitable with: Promerix, Practico, OP1650

# Accessories

## Special devices



### 95. The Beach Chair Shoulder System for shoulder arthroscopy

The Beach Chair Shoulder System gives open bilateral access for shoulder and upper arm procedures. Support belt gives extra safety for the patient on the operating side.

60525 replaces Promerix standard back rest. In Practico Beach Chair is an optional electrical back rest (factory assembled). Load capacity up to 275 kg Promerix, 225 kg Practico.

Order no.

- 1000 60525 for Promerix
- 1000 30125 for Practico (factory assembled)
- 1000 30126 Head rest for Beach Chair
- 1000 30128 Helmet type head rest
- 1000 20127 Support belt for Beach Chair

Suitable with: Promerix, Practico



### 96. Knee -Chest device for spinal procedures 60600

Positioning device with the patient in the knee - chest - position (genucubital position). The entire positioning device can be brought into the horizontal with the help of the tables motor-driven leg joint and the gas spring assisted knee section. The height of the calf section is adjustable with the help of the crank mechanism, crank can be used from both sides of the table. Support strap, buttocks - and side supports are included to the device.

Order no. 1000 60600 Knee chest device

Recommended accessories:

- Chest support D-pillow 1000 19497
- Prone head rest 1000 66111 or face pillow 1000 60751
- Trolley 1000 19285 (see the trolley on page 17)

Suitable with: Promerix



## Orthopedic extension device



### 97. Orthopedic extension device 19300

Suitable for treatment of all fractures in lower extremities. For Promerix, Practico adapter incl., for Promerix gyn. section 60425 or 60440 is needed.

Included in the delivery: Adapter to the table, dual joint extension bars (2 pcs), Support feet (2 pcs), Screw tensioners (2 pcs), Traction boots, adult model (a pair), Pelvic support, perineal post (1 pc).

Order no. 1000 19300



### 98. Short support arm 19313

With clamp 19135 for supporting healthy leg with stirrup or leg support, e.g. 20118.

Order no. 1000 19313



### 99. Leg plates 19348

For supporting legs during the anaesthesia before the final positioning.

Order no. 1000 19348




### 100. Accessory stand 19285

For 19300 and for 60600.

Order no. 1000 19285



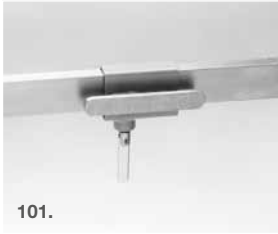

**101. Clamp for the extension bar with rail**

For wide arm 19135 for 19300.


Order no. 1000 19135



**102. Clamp for the extension bar 19134**

With 19 mm bore for 19300.


Order no. 1000 19134

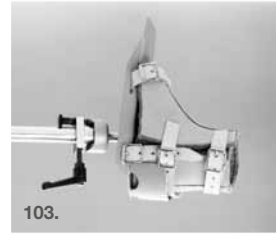

**103. Traction boots children model 19131**

Leather boots (pair), bottom plates are attached to the screw tension device with a ball joint for 19300.

Order no. 1000 19131

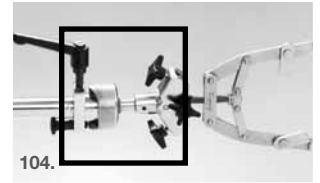

**104. Stirrup clamp with rotation 19114**

For fixing Kirchner bow (not included in the delivery) traction device. for 19300.

Order no. 1000 19114



**105. Knee crutch 19118 (m. Goepel)**

With clamp and velcro-strap, for 19300.

Order no. 1000 19118



**106. Knee support 19124**

Radiotranslucent (carbon fibre) and padded with height adjustable support, port and clamp, attaches to the extension bar, for 19300.

Order no. 1000 19124



**107. Knee support 138**

With height adjustable, telescopic support post, radiotranslucent (carbon fibre) and padded, attaches to the adapter, can be fixed on both sides of the table for 19300.


Order no. 1000 00138



**108. Leg rest 19120**

Radiotranslucent (carbon fibre) with support post and clamp, attaches to the extension bar, for 19300.


Order no. 1000 19120

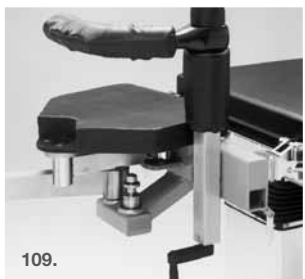

**109. Femoral counter traction post 19112**

For 19300 incl. pads, perineal posts (carbon fibre) are radiotranslucent and height adjustable.

Order no. 1000 19112

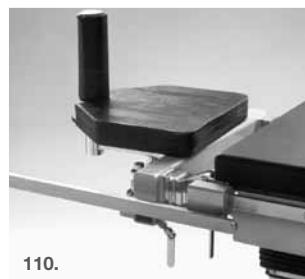


**110. Accessory rail elongation 19128**

With clamp, length 550 mm, for 19300.

Order no. 1000 19128 EUR-rail, UK-rail


1000 19128U US-rail

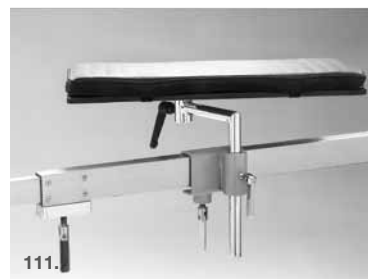

**111. Calf support 19126**

With clamp and pad, for 19300.

Order no. 1000 19126

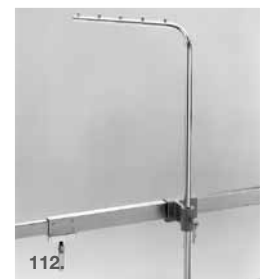

**112. Supporting post for plaster cast 19132**

For 19300.

Order no. 1000 19132

# Accessories

## Recommended accessories for the most common fields of surgery

CODE	ACCESSORIES : for EUR-rail	PROMERIX	PRACTICO	PRACTICO MANUAL	RAPIDO	OP1650
<b>ARM RESTS</b>						
100000125	Arm rest 125, with straps 110	x	x	x	x	x
100019180	Arm rest 19180, with straps 110	x	x	x	x	x
100019185	Arm Rest 19185 w.radial clamp, ball joint	x	x	x	x	x
100011801	Easy armboard 11801, with straps	x	x	x	x	x
100011802	Lateral armboard 11802	x	x	x	x	x
100010300	- for fixing : 1 pc. Radial Accessory Clamp, Universal 10300	x	x	x	x	x
<b>SUPPORTS AND STRAPS</b>						
100000107	Wrist straps 107, pair	x	x	x	x	x
100000110	Straps 110 for arm rests 125 and 19180, 19185	x	x	x	x	x
100011081	Patient restraining strap 11081	x	x	x	x	x
100010886	Leg / body strap, wide, velcro-type	x	x	x	x	x
100000120	Side support 120	x	x	x	x	x
100011651	Lateral brace 11651with round pad and clamp	x	x	x	x	x
100011652	Lateral brace 11652 with rectangular pad and clamp	x	x	x	x	x
100019129	Shoulder supports 19129, pair	x	x	x	x	x
100011295	Shoulder braces 11295, pair	x	x	x	x	x
100011043	- for fixing 2x Rail Clamp 11043	x	x	x	x	x
100020127	Support belt for Beach Chair	x	x	x	x	x
100010368	Supporting Roll, vertical	x	x	x	x	x
100010301	- for fixing Accessory Clamp, D20 mm	x	x	x	x	x
100010369	Supporting Roll, horizontal	x	x	x	x	x
100060485	Arm- / Leg Shield , incl. pad	x	x	x	x	x
<b>HEAD RESTS</b>						
100000141	Wide head rest 141	x	x	x	x	x
	-for Promerix and Practico the fixing adapter is needed 60146					
100000142	Narrow head rest 142	x	x	x	x	x
	-for Promerix and Practico the fixing adapter is needed 60146					
100000144	Neurosurgical head section 144	x	x	x	x	x
	-for Promerix and Practico the fixing adapter is needed 60146					
100018101	Head section 18101	x	x	x	x	x
	-for Promerix and Practico the fixing adapter is needed 60146					
100018143	Special head rest 18143, for ENT-surgery	x	x	x	x	x
	-for Promerix and Practico the fixing adapter is needed 60146					
100018150	Special head rest 18150, for ENT/ophtalmology	x	x	x	x	x
	-for Promerix and Practico the fixing adapter is needed 60146					
100018151	Special head rest 18151 for ophtalmology	x	x	x	x	x
	-for Promerix and Practico the fixing adapter is needed 60146					
100000151	Arm rest 151, pair	x	x	x	x	x
100030126	Special Head Rest for Beach Chair	x	x	x		
100030128	Helmet type Head Rest for Beach Chair	x	x	x		
100330010	DORO Set 1	x	x	x		
100330020	DORO Set 2	x	x	x		
100330030	DORO Set 3	x	x	x	x	x
100301100	DORO® Side Rail Adaptor, for OP1700, Practico and Promerix	x	x	x		
100060146	Adapter 60146 for special head rests in Promerix and Practico	x	x	x		
<b>FOOT RESTS</b>						
100000118	Leg rests with straps 118, pair	x	x	x	x	x
100020118	Leg rests with straps, Goebel-type, pair					
100010300	- for fixing : 2 pc. Radial Accessory Clamp, Universal 10300	x	x	x	x	x
100011181	Knee crutches 11181, with extra soft pad and straps	x	x	x	x	x
100010300	- for fixing : 2 pc. Radial Accessory Clamp, Universal 10300	x	x	x	x	x
100012551	Leg holders 12551, pair	x	x	x	x	x
100010300	- for fixing : 2 pc. Radial Accessory Clamp, Universal 10300	x	x	x	x	x
100012552	Leg support, "The Prepper" 12552	x	x	x	x	x
100010300	- for fixing : 1 pc. Radial Accessory Clamp, Universal 10300	x	x	x	x	x
100012581	Stirrup system 12581, with extra soft pad	x	x	x	x	x
100010300	- for fixing : 2 pc. Radial Accessory Clamp, Universal 10300	x	x	x	x	x
100012583	Gas spring assisted stirrups 12583, pair	x	x	x	x	x
100011043	- for fixing 2x Rail Clamp 11043	x	x	x	x	x
100010353	Arthroscopic Legholder 10353	x	x	x	x	x
100010308	- for fixing : 1 pcs. Radial Clamp 10308	x	x	x	x	x
100012530	Arthroscopic Stress Post 12530	x	x	x	x	x
100010300	- for fixing : 1 pc. Radial Accessory Clamp, Universal 10300	x	x	x	x	x
100012585	Total knee stabilizer 12585, incl. Clamp	x	x	x	x	x
100012586	Replacement pad 12586 for 12585	x	x	x	x	x
100012012	Thigh Holder (single)	x	x	x	x	x
100010200	- for fixing:1 X Accessory Clamp Universal 10200	x	x	x	x	x
100012001	Light foot rest 12001					x
100018001	Light foot rest 18001				x	
100020201	Light foot rest 20201		x	x		
100060201	Light foot rest 60201, Promerix	x	x	x		
100060482	Powered Leg Holders for Promerix	x				
<b>ANESTHESIA SCREENS AND INFUSION ROD</b>						
100000111	Anaesthesia screen 111	x	x	x	x	x
100010300	- for fixing : 1 pc. Radial Accessory Clamp, Universal 10300	x	x	x	x	x
100000116	Extension arm 116, for 111	x	x	x	x	x
100011001	Anaesthesia screen 11001, flexible	x	x	x	x	x
100010200	- for fixing:2 X Accessory Clamp Universal 10200	x	x	x	x	x
100000749	Infusion rod 749, for EUR-rail, 2 hooks	x	x	x	x	x
100000749	Infusion rod 749, for EUR-rail, 4 hooks	x	x	x	x	x
100018578	Infusion rod holder 18578, for Rapido				x	
<b>ARM AND HAND TABLES</b>						
100010390	Hand operating table 10390	x	x	x	x	x
100010392	Traction device 10392	x	x	x	x	x
100010393	Handfixation device 10393	x	x	x	x	x
100011911	Arm and hand table 11911	x	x	x	x	x
100011912	The Carpa table 11912	x	x	x	x	x

**CODE ACCESSORIES : for EUR-rail**

PROMERIX  
PRACTICO  
PRACTICO MANUAL  
RAPIDO  
OP1650

**BOWLS**

10000158	Bowl 158, 4 L + attachment, fixed under the base						x
100181581	Bowl 181581, 6 L , fixed under the seat section without drainage					x	
100181582	Bowl 181582, 6 L , fixed under the seat section with drainage					x	
100201501	Bowl 201501, 6 L, without drainage	x	x	x	x		
100201502	Bowl 201502, 6 L, with drainage	x	x	x	x		
100201571	Bowl 201571 6 L + attachment, fixed to the side without drainage	x	x	x	x		x
100201572	Bowl 201572 6 L + attachment, fixed to the side with drainage	x	x	x	x		x
100061580	Holder for 6 l Bowl, , fixed under the base, Promerix	x					
100061581	Bowl 6 l without Drainage + Holder, fixed under the base, Promerix	x					
100061582	Bowl 6 l with Drainage + Holder, fixed under the base, Promerix	x					
100018159	Screen 18159 for bowls, 6 L	x	x	x	x		x

**OTHER ACCESSORIES**

10000126	Cassette tray 126	x	x	x			x
100018126	Casette tray 18126					x	
100018153	Gynecologic extension 150 mm					x	
100018168	Push handles to foot end	x	x	x	x		x
100018173	Side rails, Chrome plated, autolock	x	x	x	x		x
100018174	Side rails, epoxy coated, autolock	x	x	x	x		x
100018183	Wire basket for utilities					x	
100018189	Fixation parts for oxygen bottle 5 L					x	
100018572	Paper roll holder, for Rapido					x	
100020593	PUSH BAR, head end	x	x	x			
100060850	Foot Control device (height, trend, tilt)	x	x	x			
100020114	Foot Rest / Table extension, without the mattress	x	x	x	x		x
100020114	Foot Rest / Table extension, with IS-molded mattress with pin screws	x	x	x	x		x
100020114	Foot Rest / Table extension, with VEF-welded visco mattress with velcro	x	x	x	x		x
100010659	Table width extender 10659, 500 mm	x	x	x			
100060484	Divided leg plates, 4-sectional	x					
100060525	Beach Chair Back Rest for shoulder athroscopy	x					
100019280	Back and kidney elevator 19280 with mattress	x	x	x			x

**CLAMPS**

100010200	Accessory Clamp 10200, Universal D16	x	x	x	x		x
100010300	Radial Accessory Clamp 10300, Universal D16-D20	x	x	x	x		x
100010308	Radial Clamp	x	x	x	x		x
100011042	Easylock Socket 11042, max. D 18 mm	x	x	x	x		x
100011043	Rail clamp 11043	x	x	x	x		x
100019128	Accessory rail extension 19128	x	x	x	x		x
100010355	Accessory Extension Rail, L=400 mm	x	x	x	x		x
100010356	Accessory rail extension, angeled	x	x	x	x		x
100010359	Angular accessory rail, L=250 mm	x	x	x	x		x

**PADS AND MATTRESSES**

100019400	Transfer Mattress	x	x	x	x		x
100019497	D-Pillow, visco-elastic, includes fixing straps,	x	x	x	x		x
100060750	Mattress set for prone position	x	x	x	x		x
100060751	Prone Head pillow	x	x	x	x		x
100040201	Donut head pad 40201	x	x	x	x		x
100040402	Gel pad, small 40402	x	x	x	x		x
100040204	Horseshoe head pad 40204	x	x	x	x		x
100040308	Large armboard pad 40308, size: 61x16,5x1 cm	x	x	x	x		x
100040617	Patient positioner, flat bottom 40617	x	x	x	x		x
100040100	O.R. Overlay, medium, size: 117x50x1,3 cm	x	x	x	x		x
100040105	O.R. Overlay, small, size: 50x50x1,3 cm	x	x	x	x		x
100040202	Donut head pad 40202, Pediatric	x	x	x	x		x
100040205	Horseshoe head pad 40205, Pediatric	x	x	x	x		x
100040213	Contoured Head Pad, size: 18x18x8 cm	x	x	x	x		x
100040300	Armboard Pad, short, size: 38x10x0,6 cm	x	x	x	x		x
100040301	Armboard Pad, long, size: 61x11,5x1 cm	x	x	x	x		x
100040601	Axillary Roll Cover, with hook&loop, size: 43x30,5x1,3 cm	x	x	x	x		x
100066002	Closed head ring / 140 x 60 x 35mm/	x	x	x	x		x
100066003	Closed head ring / 200 x 75 x 45mm/	x	x	x	x		x
100066012	Open head ring / 140 x 60 x 35mm/	x	x	x	x		x
100066013	Open head ring / 200 x 75 x 45mm/	x	x	x	x		x
100066021	Table pad / 500 x 500 x 10mm	x	x	x	x		x
100066022	Table pad / 1170 x 500 x 10mm	x	x	x	x		x
100066096	Patient positioner / 380 x 125 x 63mm	x	x	x	x		x
100066042	Chest roll with flat bottom surface / 405x150x100 mm	x	x	x	x		x
100066072	Armboard pad /600 x 125 x 20mm	x	x	x	x		x
100066031	Heel support pads, pair / 200 x 110 x 70mm	x	x	x	x		x
100066801	Crutch stirrup pads-pair / 315x220x10 mm	x	x	x	x		x
100066081	Contoured headrest / 190x190x60 mm	x	x	x	x		x
100066110	Head pad / 270x230x70 mm	x	x	x	x		x
100066111	Prone headrest / 280 x 230 x 145mm	x	x	x	x		x

**SPECIAL DEVICES**

100000154	Proctology attachment 154						x
100030154	Proctology attachment for Promerix, Practico	x	x	x			
100019220	Accessory stand 19220 (incl. 4 wire baskets)	x	x	x	x		x
100060500	Carbon fibre backrest, long, 950 mm	x	x	x			
100060501	Mattress for long carbon fiber section, 950 mm, VEF 80 mm	x	x	x			
100060510	Arm rest, Carbon fiber	x	x	x			
100060511	Pad for carbon fiber arm rest, VEF 80 mm	x	x	x			
100060600	Knee Chest device for spinal surgery, Promerix	x					
100019300	Orthopaedic Extension Device, Promerix / Practico	x	x	x			
	Extension bars, 2 pcs.						
	Support feet, 2 pcs.						
	Screw tensioners, 2 pcs.						
	Traction boots (adults'), 2 pcs.						
	Pelvic support + perineal post, 1 + 1 pc.						

**OPTIONAL ACCESSORIES FOR 19300**

100000138	90° knee support 138	x	x	x			
100019112	Femoral counter traction post 19112	x	x	x			
100019114	Stirrup clamp with rotation 19114	x	x	x			
100019118	Knee crutch (model Göpel) 19118	x	x	x			
100019120	Leg rest with padding 19120	x	x	x			
100019124	Knee support 19124	x	x	x			
100019126	Calf support 19126	x	x	x			
100019128	Accessory rail extension 19128	x	x	x			
100019131	Traction boots, children's 19131	x	x	x			
100019132	Supporting post for plaster cast 19132	x	x	x			
100019134	Clamp 19134 for the extension bar (Ø 19 mm hole)	x	x	x			
100019135	Rail clamp 19135 for wide arm	x	x	x			
100019285	Cart for Orthopaedic device and Knee Chest device	x	x	x			
100019313	Short support arm for 19300	x	x	x			
100019348	Leg plates for Othopaedic device 19300	x	x	x			


# Solutions for efficient physical patient flow and well-being

Merivaara provides a wide range of hospital-grade furniture such as operating tables, medical lights, examination tables, trolleys and stretchers for transportation and day surgery, as well as birthing and patient beds. Additionally, our new concepts include both revolutionary integrated management systems and after-sales services.

Our versatile solutions and product offering enable convenient and cost-effective physical patient flow within the hospital, in day surgery and health clinics, as well as in nursing homes and home care. Our priority is to enhance comfort and well-being of both patients and their caregivers.

Established in 1901, Merivaara has more than 100 years of experience in designing and manufacturing hospital-grade furniture. Today our products are highly appreciated by users in more than 120 countries thanks to their ease-of-use, durability and ergonomic design.

Merivaara's medical products are CE-labelled and the company's quality management system complies with EU directives for medical device. It is certified by DNV according to ISO 9001:2008 and ISO 13485:2003 standards. The company's environmental management system is ISO 14001:2004 certified.

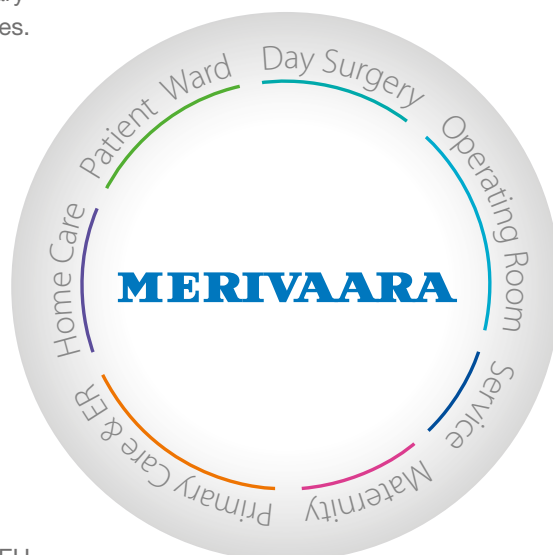

**Merivaara Corp.**  
**Puustellintie 2, FI-15150 Lahti, FINLAND**  
**Tel. +358 3 3394 611, fax +358 3 3394 6144**  
**merivaara@merivaara.com**  
**www.merivaara.com**

The status of the information, specifications and illustrations in this publication is indicated by the date of publication. Whilst the greatest care has been taken over its compilation, Merivaara Corp. is not responsible for any errors or omissions. Merivaara Corp. reserves the right to make changes to the technology, features, specifications and design of the equipment and models without notice. Please check the Merivaara website for the latest updated specifications and other information.

All trademarks are Merivaara property unless otherwise stated. Patents granted or pending.